

waterpik®
waterflosser™

Congratulations!

You are now on your way
to better oral health!

✓ Register and Save

Register your product at *waterpik.com/registration to receive over \$50 in savings (includes 50% off your next order of assorted tips + tip storage case)

✓ Write a Review

Creating innovative Waterpik® products that inspire better oral health is our passion. We invite you to *leave a review by sharing your thoughts on your new Water Flosser on waterpik.com or the retailer's website.

*US customers only

✓ We're Here to Help

We hope you are loving your new Water Flosser. If for any reason you are not satisfied or have any product questions, please contact our fast and friendly Customer Service Team based in Fort Collins, Colorado. your new Water Flosser is backed by a Two Year manufacturer's warranty.

Call: 1-800-525-2774 (USA)
Monday-Friday
7 am - 5 pm MST

Chat or E-mail
service@waterpik.com



**Quick Start
Guide Inside!**



Register and Save!

Save over \$50!

Register at *Waterpik.com/Registration to receive exclusive discount code.

Offer includes
50% off
your next order
of assorted
tips + tip
storage case.



We'd Love to Hear from You!



Connect with Us
@Waterpik



Call: 1-800-525-2774 (USA)
www.waterpik.com

Waterpik® 14 Day Challenge App



Make water flossing a habit! Try our new mobile App for the best water flossing experience.




Download on the
App Store

GET IT ON
Google play

Quick Start Guide


GETTING STARTED

- 1 Turn power switch **OFF(0)**, plug power cord into outlet, and insert the charging cord into the unit.

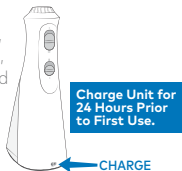


- 2 Insert tip firmly until it clicks into place.

Caution:
Do not press tip eject button while unit is running.

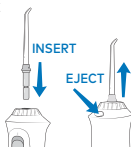

- 3 Unplug charging cord from unit prior to filling with water and during use. Lift the flip top on reservoir and fill with lukewarm water. Refill as needed.

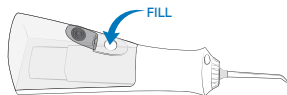

WATER FLOSS

- 4 For first time use, adjust the pressure setting to the **LOW(I)** setting.



- 5 Lean over sink and close lips enough to prevent splashing, while still allowing water to flow from mouth into the sink.



Turn unit **ON(I)**.

TECHNIQUE

- 6 Aim water at the gumline at a 90 degree angle.

Follow the gumline and pause briefly between teeth.



NEED HELP?

Visit Waterpik.com/WelcomeWF
To watch a How To Video and find additional product support.

**Call Us! 1-800-525-2774 (USA)
1-888-226-3042 (CAN)**

TO AVOID MESS: KEEP TIP IN MOUTH WHILE UNIT IS ON.