
**waterpik®**
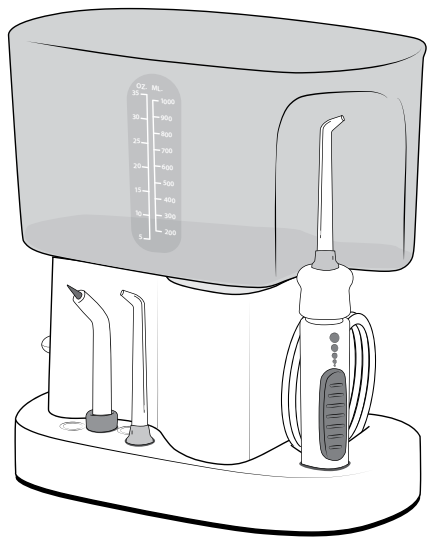
## **Waterpik® Water Flosser Model WP-70 Series**

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## IMPORTANT SAFEGUARDS

When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

### READ ALL INSTRUCTIONS BEFORE USING.

## DANGER:

### To reduce the risk of electrocution:

- Always unplug product after using.
- Do not handle plug with wet hands.
- Do not place in or drop into water or other liquid.
- Do not use while bathing.
- Do not place or store product where it can fall or be pulled into a tub, or sink.
- Do not reach for a product that has fallen into water. Unplug immediately.
- Check the power supply cord for damage before the first use and during the life of the product.

## WARNING:

### To reduce the risk of burns, electrocution, fire, or injury to persons:


- Do not plug this device into a voltage system that is different from the voltage system specified on the device.
- Do not use this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into any water. Contact Waterpik International, Inc. or your local distributor which can be found at [www.waterpik.com/intl](http://www.waterpik.com/intl).
- Do not direct water under the tongue, into the ear, nose or other delicate areas. This product is capable of producing pressures that may cause serious damage in these areas.
- Use this product only as indicated in these instructions or as recommended by your dental professional.
- Only use tips and accessories recommended by Water Pik, Inc.

- Do not drop or insert any foreign object into any opening or hose.
  - Keep electrical cord away from heated surfaces.
  - Do not use outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.
  - Do not use iodine, saline (salt solution), or water insoluble concentrated essential oils in this product. Use of these can reduce product performance and will shorten the life of the product. Use of these products may void warranty.
  - Remove any oral jewelry prior to use of this product.
  - Do not use if you have an open wound on your tongue or in your mouth.
  - If your physician or cardiologist has advised you to receive antibiotic premedication before dental procedures, you should consult your dentist before using this instrument or any other oral hygiene aid.
  - If the Pik Pocket™ tip (not included in all models) separates from the shaft for any reason, discard tip and shaft and replace with new Pik Pocket™ tip.
  - Do not use this product for more than 5 minutes in each two-hour period.
- Closely instruct and supervise children 8 years and above and individuals with special needs in the proper use of this product.
  - Children should be supervised to ensure that they do not play with the appliance.
  - There are no consumer serviceable parts in this appliance and it does not require maintenance.
  - Fill reservoir with warm water or other dental professional recommended solution only.


### **SAVE THESE INSTRUCTIONS**

Product Description	5
Getting Started	6
Cleaning and Troubleshooting Your Waterpik® Water Flosser	8
Limited Two-Year Warranty	9

1. Reservoir
2. Pump Inlet
3. Base
4. ON/OFF
5. Tip
6. Rotating direction control knob
7. Handle
8. Pressure Control Switch
9. Cord Storage (under base)
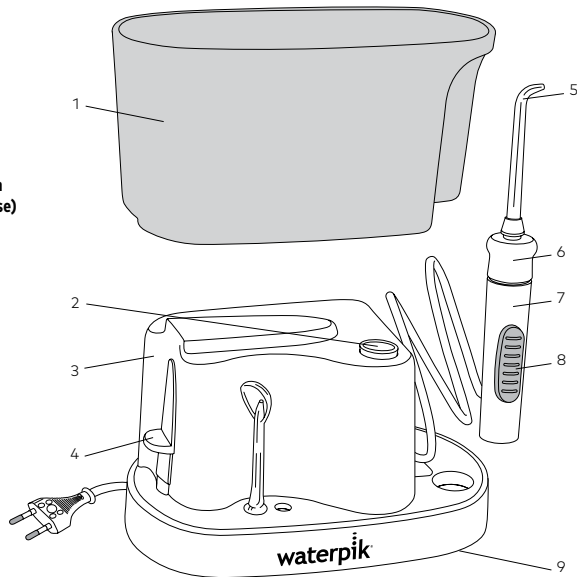
#### NOT ALL TIPS INCLUDED IN ALL MODELS



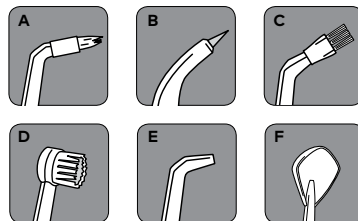
\*Replace Every 3 months \*\*Replace Every 6 months

#### TIPS

- A. Orthodontic Tip\***  
- Braces/General Use
- B. Pik Pocket™ Tip\***  
- Periodontal Pockets/Furcations
- C. Plaque Seeker™ Tip\***  
- Implants/Crowns/Bridges/  
Retainers/General Use
- D. Toothbrush Tip\***  
- General Use
- E. Classic Jet Tip\*\***  
- General Use
- F. Tongue Cleaner\*\***  
- Fresher Breath

Replacement tips can be purchased from your local distributor, which can be found at [www.waterpik.com/intl](http://www.waterpik.com/intl).
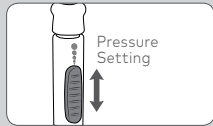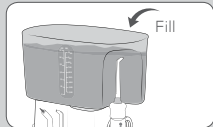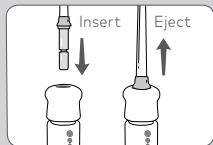
### Powering the Unit

Plug the cord into an electrical outlet. If the outlet is controlled by a wall switch, be sure it is turned on.

### Inserting and Removing Tips

Pull the knob on the top of the handle toward the hose and insert tip into the center of the knob. The colored ring will be flush with the end of the knob if the tip is correctly locked into place.

To remove tip from the handle, pull the knob on the top of the handle toward the hose and pull tip from handle.

### Preparing the Reservoir

Remove reservoir from base. Fill the reservoir with lukewarm water. Replace the reservoir over the base and press down firmly.

### Prime Unit (First Time Use Only)

Set pressure control dial on base to high, point the handle and tip into sink. Turn unit ON (I) until water flows from the tip. Turn unit OFF (O).

### Adjusting the Pressure Setting

Slide the pressure control switch located on the handle to the lowest setting (away from the tip as far as possible) for first-time use. Gradually increase pressure over time to the setting you prefer – or as instructed by your dental professional.

### Recommended Technique

For best results, start with the back teeth and work toward the front teeth. Glide tip along gumline and pause briefly between teeth. Continue until you have cleaned the inside and outside of both the upper and lower teeth.

Lean over sink. Direct the stream at a 90-degree angle to your gumline with tip in mouth. **Slightly close lips to avoid splashing**, but allow water to flow freely from mouth. Set pressure at desired setting and turn the unit ON (I). Begin cleaning.

### Pause Control

You can temporarily stop the flow any time by pressing the pause button on the handle.

## When Finished

Turn the unit off (o). Empty any liquid left in the reservoir. Unplug the unit.

**Note:** if you used a mouthwash or other solution besides water, see "using Mouthwash or other solutions" section.

## NOTE: NOT ALL TIPS INCLUDED WITH ALL MODELS.

### TIP USAGE

#### Pik Pocket™ Tip

The Pik Pocket™ tip is specifically designed to deliver water or anti-bacterial solutions deep into periodontal pockets.

To use the Pik Pocket™ tip, **set the unit to the lowest pressure setting**. Place the soft tip against a tooth at a 45-degree angle and gently place the tip under the gumline, into the pocket.

#### Tongue Cleaner Tip

Set to lowest pressure setting, place tip in the center/middle of your tongue about half way back. Pull forward with light pressure. Increase pressure as you prefer.

#### Orthodontic Tip


To use the orthodontic tip, gently glide tip along gumline, pausing briefly to lightly brush area between teeth and all around orthodontic bracket, before proceeding to the next tooth.

#### Toothbrush Tip

Place the toothbrush tip in mouth with the brush head on the tooth at the gumline. The toothbrush tip can be used with or without toothpaste. Using a **light pressure** (bristles should not bend) massage the brush back and forth with very short strokes – much like you would with a manual toothbrush.

#### Plaque Seeker™ Tip

To use the Plaque Seeker™ tip, place the tip close to the teeth so that the bristles are gently touching the teeth. Gently glide tip along gumline, pausing briefly between teeth to gently brush and allow the water to flow between the teeth.
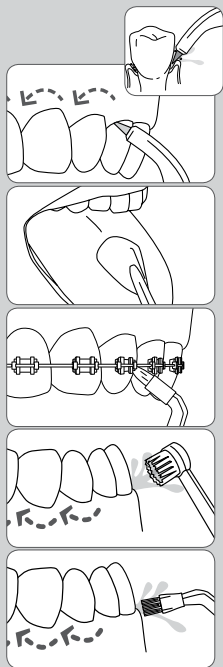
### Using Mouthwash and Other Solutions

Your Waterpik® Water Flosser can be used to deliver mouthwash and antibacterial solutions. After using any special solution, rinse unit to prevent clogging by partially filling reservoir with warm water and running unit with tip pointed into sink until unit is empty.

### Cleaning

Before cleaning, unplug from the electrical outlet. Clean the product when necessary by using a soft cloth and mild non-abrasive cleanser. The water reservoir is removable for easy cleaning and is top rack dishwasher safe.

**NOTE:** Remove reservoir valve before placing reservoir in dishwasher.



### Removing Hard Water Deposits

Add 2 tablespoons of white vinegar to a full reservoir of warm water. Set pressure control dial to high, point the handle and tip into sink. Turn unit ON (I) and run until reservoir is empty. Rinse by repeating with a full reservoir of clean warm water.

### Leakage

If leaking due to a damaged hose, replacement hoses and handles can be purchased from your local distributor, which can be found at [www.waterpik.com/intl](http://www.waterpik.com/intl).

If leakage occurs from the reservoir, ensure the valve is right side up or contact your local distributor for a replacement valve.
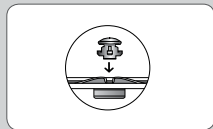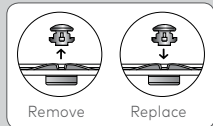
### **Service Maintenance**

Waterpik® Water Flossers have no consumer-serviceable electrical items and do not require routine service maintenance. Some replacement parts and accessories for the Water Flosser are also available from your local distributor which can be found at [www.waterpik.com/intl](http://www.waterpik.com/intl).

Refer to the serial and model numbers in all correspondence. These numbers are located on the bottom of your unit. Still have questions? Go to [www.waterpik.com/intl](http://www.waterpik.com/intl) to find your local distributor.

Water Pik, Inc. warrants to the original purchaser/owner of this product that it is free from defects in materials and workmanship for two years from date of purchase. Save your receipt as your proof of purchase date. Water Pik, Inc. reserves the right to request the product to be returned for analysis to confirm the claim falls within our warranty policy. We will replace any part of the product, which in our opinion is defective, provided the product has not been abused, misused, altered or damaged after purchase and was used according to these instructions. Repair or replacement services provided for by the guarantee will be carried out by the local distributor, whose details may be found on [[www.waterpik.com/intl](http://www.waterpik.com/intl)]. Upon request of the consumer, a copy of the guarantee complete with the local distributor's information will be provided by email or mail. This limited warranty excludes accessories or consumable parts such as tips, etc.

This warranty gives you specific legal rights. You may also have other rights that vary depending on the law in your place of residence. To locate your local distributor, go to [www.waterpik.com/intl](http://www.waterpik.com/intl).

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©2019 Water Pik, Inc.
**Water Pik, Inc., a subsidiary
of Church & Dwight Co., Inc.**
SOFIBEL SAS
Church & Dwight UK LTD.,